

Compact, automated nucleic acid extraction system

SEEPREP³²[™]

SEEPREP³²[™], the automated extraction system, can simultaneously extract up to 32 specimens within just 30 minutes, which guarantees to enhance productivity and efficiency of your laboratory.

○ Application of a wide range of specimen types

- Simultaneous extraction of nucleic acid (DNA/RNA) from various types of specimens

- Aspirate (Nasopharyngeal)	- Plasma	- Stool (Cary-Blair)
- BAL (Bronchoalveolar lavage)	- Saliva	- Swabs (Nasopharyngeal, Vaginal, Cervical, Urethral, Rectal)
- CSF	- Serum	- Urine
- LBC (Liquid based cytology)	- Sputum	- Whole blood

○ Short turnaround time

- Extraction of 1-32 specimens in less than 30 minutes

○ Superior user convenience

- A simple interface and a remote tablet operation
- Ready-to-use reagent cartridge
- No loss of reagents when using any of 1 to 32 specimens per run



(STARMag 96 ProPrep / STARMag 96 ProPrep plus)



SEEPREP³²[™]

Compact, automated nucleic acid extraction system

More efficient and effective extraction system

1 Easy operation

- Easy to monitor the whole process
- Operates by tablet as a remote control mode
- Allows the operator to create a new protocol by tablet

2 Qualified performance


- Highly purified extracted nucleic acid
- Magnetic bead transfer-based extraction

3 Ready-to-use reagent cartridge

- Prefilled extraction reagents for user convenience
- Only one protocol covering various types of specimens

4 Prevention of contamination


- Internal sterilization with built-in UV light

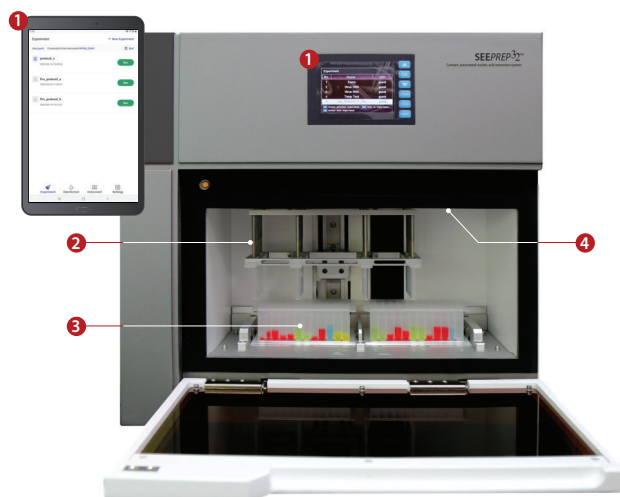

The principle of the magnetic bead transfer-based nucleic acid extraction

SEEPREP³²[™] takes advantages of the magnetic rods that automatically transfer the magnetic beads-nucleic acid mixture from one well to another, then process the cell lysis, nucleic acid adsorption, washing and elution.

- Minimizing PCR inhibitors
- Rapid extraction
- Reducing consumables

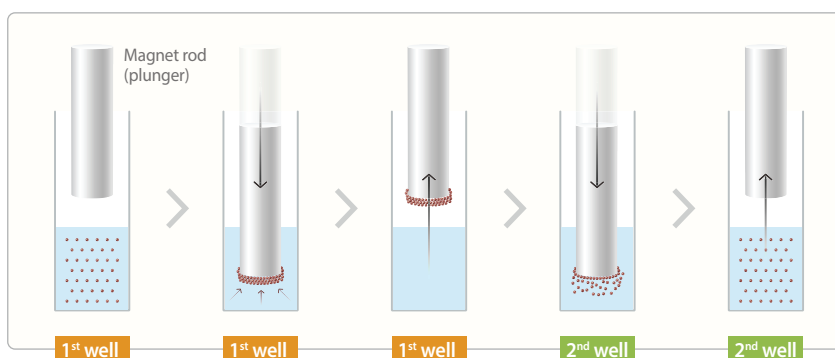

Technical specification

Extraction system method	Magnetic bead transfer
Processing capacity	1 to 32 specimens per run
Turnaround time	Within 30 minutes for 32 specimens
Hands-on time	5-10 minutes
Convenience function	Wireless control tablet
Power input	100-240VAC, 50-60Hz
Power consumption	600W

Dimension (W x D x H)	440 x 430 x 445 mm
Weight	31.5kg
Operating temperature	20-30°C
Relative humidity	10% - 80% R.H. (non-condensing)
Magnetic rod	32 ea
Sample volume	200 µL
Elution volume	100 µL

Ordering information

Category	Cat. No.
Instrument SEEPREP³²[™]	SG71100
Extraction reagent STARMag 96 ProPrep	EX00009P / EX00009T
Extraction reagent STARMag 96 ProPrep C	EX00017P / EX00017T
Extraction reagent STARMag 96 ProPrep plus*	EX00010P / EX00010T

* only for Whole Blood, Serum, Plasma



Seoul, Republic of Korea / Tel : (+82) 2 2240 4000 / E-mail : info@seegene.com

seegene.com

BRAZIL
Contagem, MG, Brazil
(+55) 31 25153003
contato@seegenebrazil.com.br

CANADA
Toronto, Canada
(+1) 800 964 5680
canada@seegene.com

GERMANY
Düsseldorf, Germany
(+49) 211 83831 600
sgg@seegene.com

ITALY
Genova, Italy
(+39) 010 667 1796
info@arrowdiagnostics.it

MEXICO
México City, México
(+52) 55 5132 6490
marketing.mx@seegene.com

MIDDLE EAST
Dubai, UAE
(+971) 4 558 7110
sgme@seegene.com

USA
California, USA
(+1) 925 448 8172
info.seegeneusa@seegene.com