


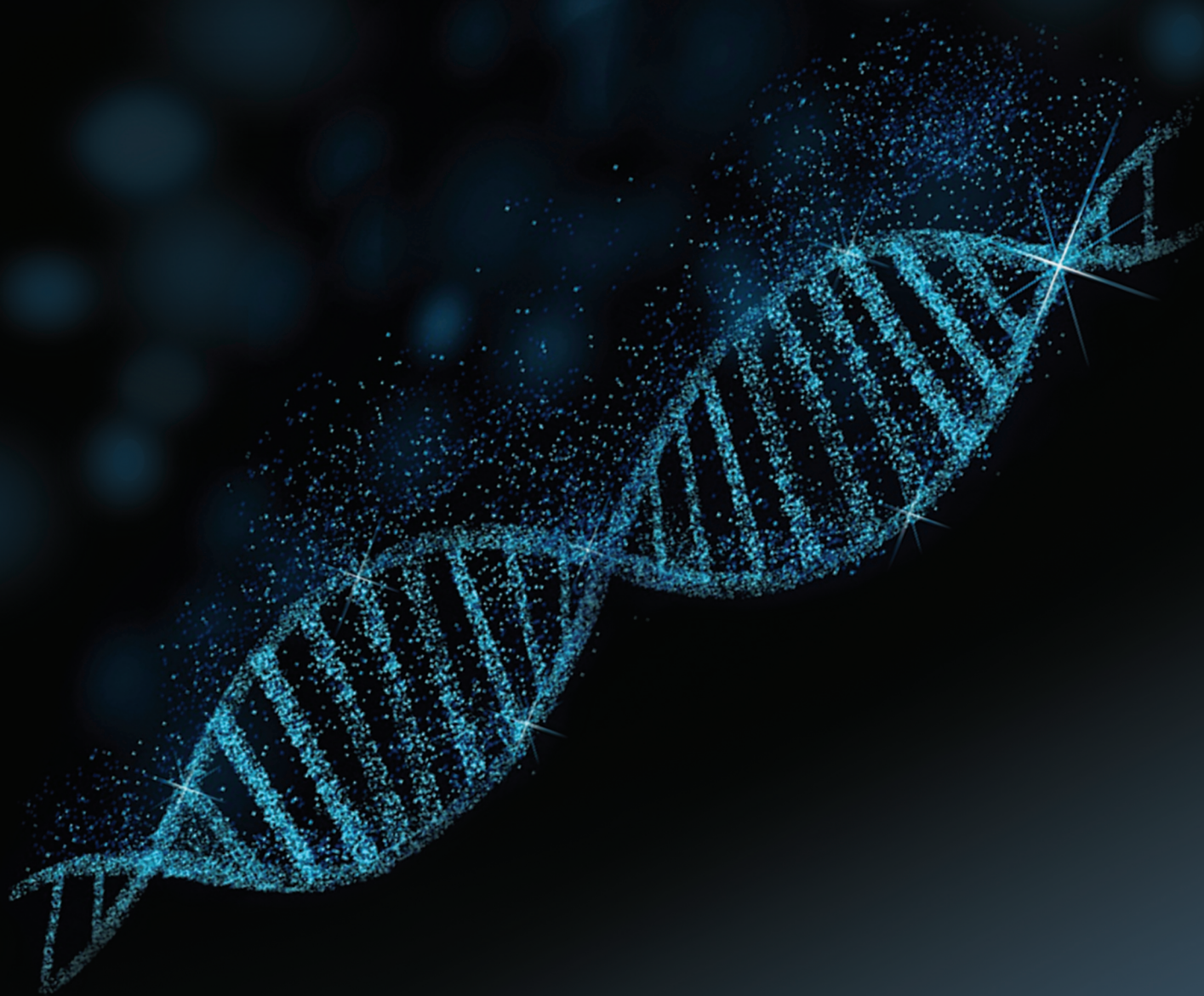


PRODUCTS FOR
MOLECULAR ONCOLOGY LABS

WHO WE ARE

EntroGen offers end-to-end solutions in molecular diagnostics for oncology and hematology by providing hardware, reagents and software. Starting from automated nucleic acid extraction and culminating in genomic interpretation with our proprietary software embedded for key decision support, our solutions are designed to be accurate, economic and scalable. The company mission is to develop diagnostic solutions for therapeutics and help drive the advancement of targeted therapies in oncology.

Through ongoing research and innovation, we are gaining worldwide recognition as leaders in the area of molecular diagnostics providing accurate and reliable products that our clients can depend on. Our aim is to provide laboratory professionals with the tools and support they need to enhance their services while managing time and cost.

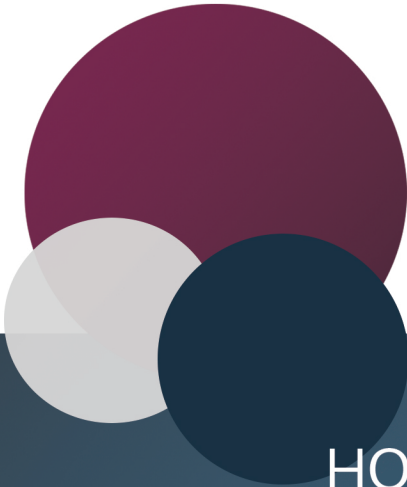

HOW WE CAN SUPPORT YOU



Internal Review

Review of current molecular diagnostic services & procedures.



Evaluation

Identify areas to improve time, cost & efficiency.



Implementation

Reagents are compatible with most leading qPCR/NGS instruments.



Automation

Fully automated sample prep and plate setup.



Software

Complimentary software for automated analysis and reporting.



Support

Prompt, knowledgeable and reliable technical support.




REAL-TIME PCR

Real-time polymerase chain reaction (rtPCR) is a well-established, rapid and highly sensitive molecular diagnostic technique used to detect specific gene alterations in DNA and RNA.

Using specific primers and probes, EntroGen rtPCR assays amplify specific targeted mutations present in a sample nucleic acid sequence. They can be easily implemented in any cancer diagnostics laboratory for detection of clinically relevant SNVs, small indels and/or gene fusions.

FEATURES

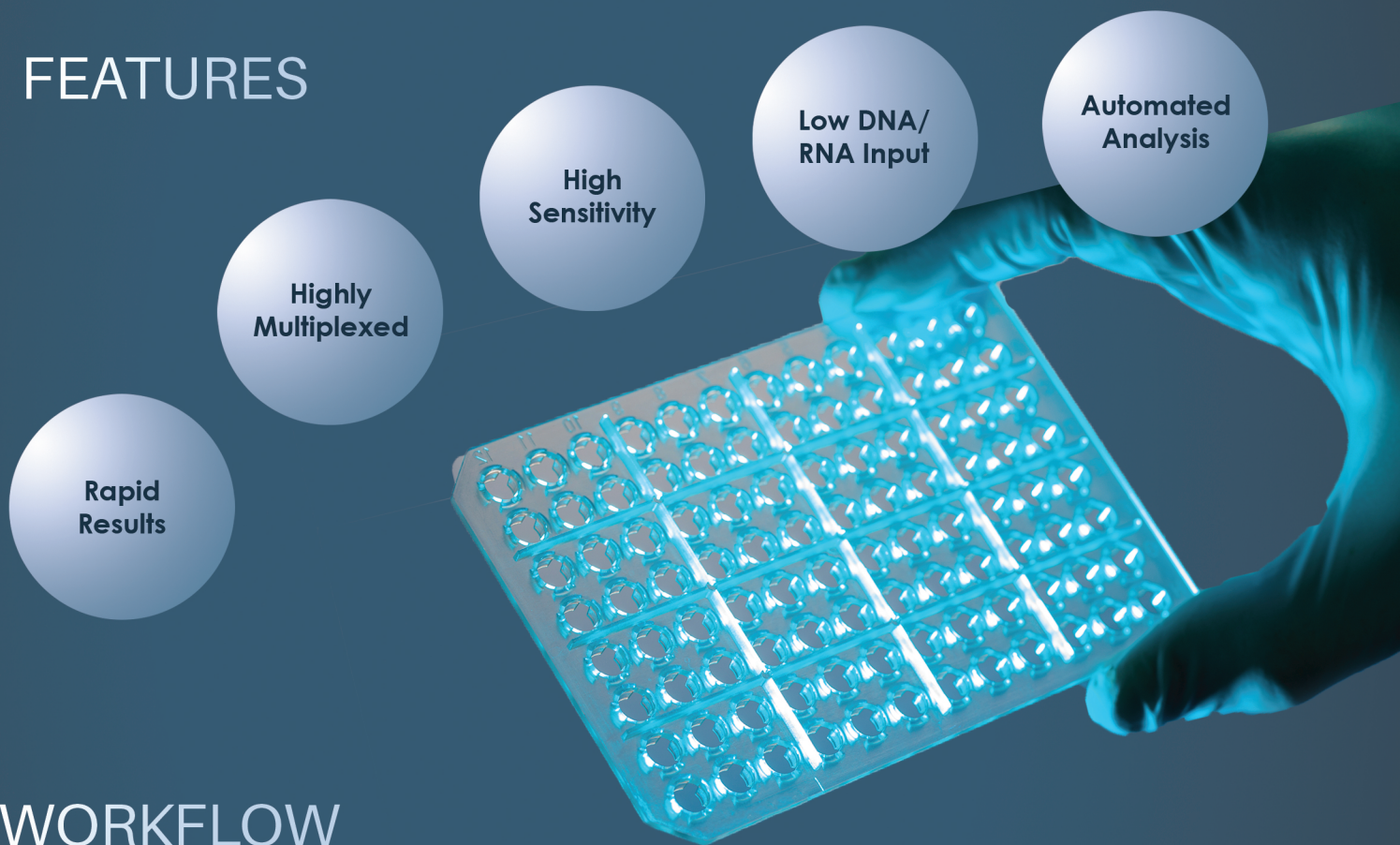

WORKFLOW



CRCdx® RAS Mutation Detection Kit

Companion Diagnostic for Vectibix®


IVD

Features

- Simple Workflow - 90 minutes from PCR setup to results
- Multiplexed - detects 35 mutations in 8 reactions
- Highly sensitive and specific - more sensitive than NGS
- Automated interpretation - includes analysis software



The most comprehensive real-time PCR companion diagnostic for RAS mutation testing

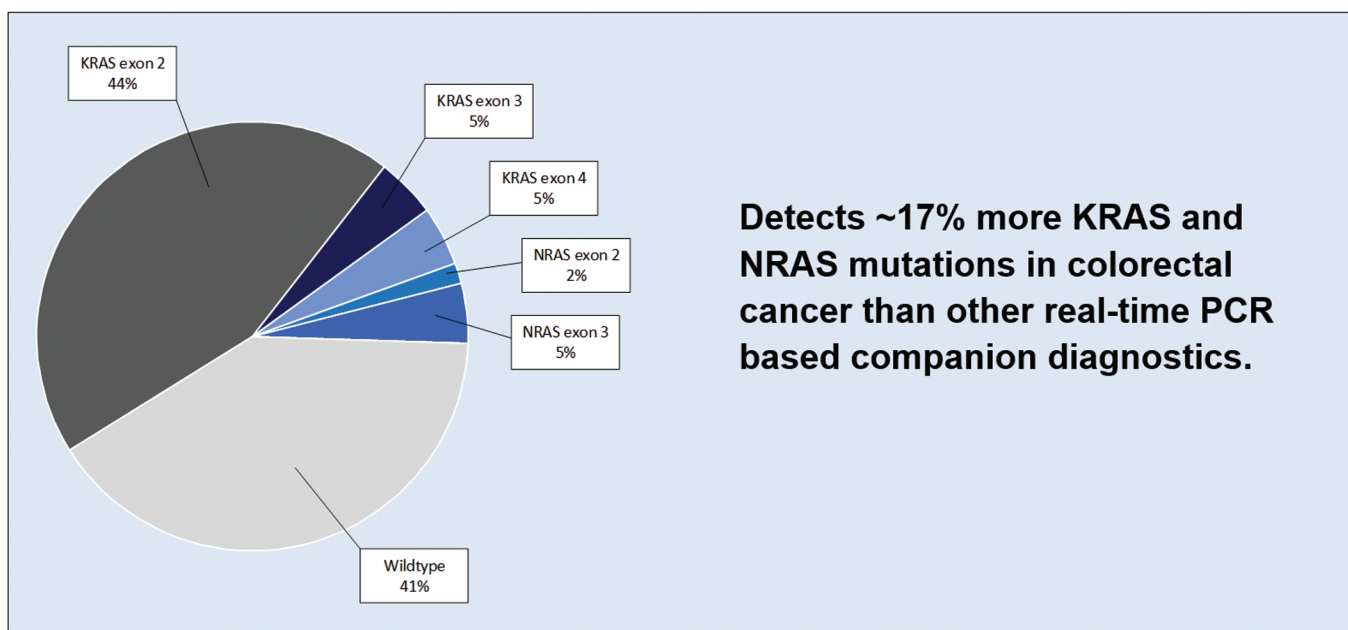

Clinical Performance

- 99.6% PPA
- 92.0% NPA
- 96.2% OPA



SOLID TUMOR ASSAYS

Colorectal Cancer

COLORECTAL CANCER MUTATION DETECTION PANEL

Targets	KRAS/NRAS G12/13, Q61, A59, K117, A146 PIK3CA E542/545, H1047 AKT1 E17K & BRAF V600	Sensitivity	As low as 1.0%
---------	---	-------------	----------------

KRAS/BRAF MUTATION ANALYSIS KIT

Targets	KRAS G12/13, Q61, K117, A146 BRAF V600E	Sensitivity	As low as 0.25%
---------	--	-------------	-----------------

RASX Mutation Analysis Extension Kit available with this assay

PIK3CA MUTATION ANALYSIS KIT

Targets	C420,E542,E545, Q546,H1047	Sensitivity	As low as 1.2%
---------	----------------------------	-------------	----------------

MLH1 METHYLATION DETECTION KIT

Targets	Methylation detection of MLH1 promoter regions C and D	Sensitivity	6.25% Methylation
---------	--	-------------	-------------------

MICROSATELLITE INSTABILITY DETECTION KIT

Targets	KIT,MSH2,SLC7A8, ZNF-2, MAP4K3	Sensitivity	As low as 5%
---------	--------------------------------	-------------	--------------

RAS MUTATION ANALYSIS KIT*

Targets	KRAS/NRAS G12/13, Q61, A59, K117, A146	Sensitivity	As low as 0.32%
---------	--	-------------	-----------------

*An add-on reaction that distinguishes KRAS G12C from other codon 12 mutations is available upon request.

KRAS MUTATION ANALYSIS KIT

Targets	G12/13, Q61, K117, A146	Sensitivity	As low as 0.25%
---------	-------------------------	-------------	-----------------

RASX Mutation Analysis Extension Kit available with this assay

NRAS MUTATION ANALYSIS KIT

Targets	G12/13, Q61, A146T	Sensitivity	As low as 0.25%
---------	--------------------	-------------	-----------------

RASX Mutation Analysis Extension Kit available with this assay

RASX MUTATION ANALYSIS EXTENSION KIT

Add-on kit for KRAS and NRAS Mutation Analysis Kits

Targets	KRAS A59 T/E/G , NRAS A59 T/D & K117R	Sensitivity	As low as 2.5%
---------	---------------------------------------	-------------	----------------



Lung Cancer

EGFR MUTATION ANALYSIS KIT

Targets	T790M, Ex19Del, L858R, L861Q, S768I, G719X, Ex20ins	Sensitivity	As low as 0.25%
---------	---	-------------	-----------------

EGFRX Mutation Analysis Extension Kit available with this assay

EGFRX MUTATION ANALYSIS EXTENSION KIT

Add-on kit for EGFR Mutation Analysis Kit

Targets	L858M, L861R & E709 A/D/G	Sensitivity	As low as 0.1%
---------	---------------------------	-------------	----------------

LUNG CANCER RNA PANEL

Targets	ALK, ROS1, RET fusions MET exon 14 skipping mutation	Sensitivity	As low as 10 copies
---------	---	-------------	---------------------

NTRK RNA FUSION GENE DETECTION KIT

Targets	NTRK 1,2,3.	Sensitivity	As low as 25 copies
---------	-------------	-------------	---------------------



Glioblastoma


MGMT METHYLATION DETECTION KIT*

Targets	MGMT promoter methylation	Sensitivity	3.125% Methylation
---------	---------------------------	-------------	--------------------

*Bisulfite Conversion Kit included

IDH1/2 MUTATION DETECTION KIT


Targets	IDH1 R100, R132 IDH2 R140, R172	Sensitivity	As low as 1%
---------	------------------------------------	-------------	--------------

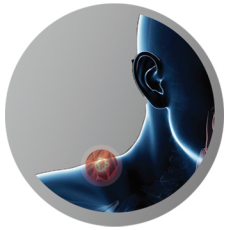

Melanoma

BRAF CODON 600 MUTATION ANALYSIS KIT II			
Targets	BRAF V600 (E, E2, K, R, D, M, G)	Sensitivity	As low as 0.25%

c-KIT MUTATION DETECTION KIT			
Targets	Exons 9, 11, 13, 17	Sensitivity	As low as 1%

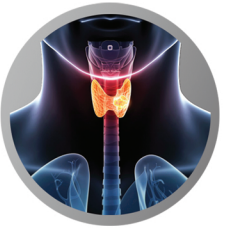

Thyroid Cancer


THYROID CANCER MUTATION DETECTION KIT			
Targets	BRAF V600E KRAS G12/13 NRAS Q61 HRAS G12/13, Q61R	Sensitivity	As low as 0.25%

THYROID CANCER FUSION GENE DETECTION KIT			
Targets	RET/PTC1 RET/PTC3 PAX8/PPAR γ	Sensitivity	As low as 0.78%



Gastrointestinal Stromal Tumors (GIST)

GIST MUTATION DETECTION KIT			
Targets	cKIT exons 9, 11, 13, 17 PDGFRA exon 18	Sensitivity	As low as 1%

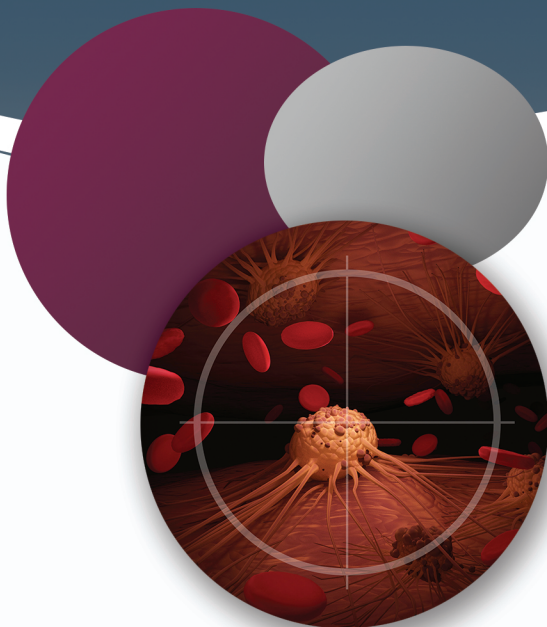

LIQUID BIOPSY ASSAYS

The utility of cell-free tumor DNA (ctDNA) from plasma has been demonstrated as an alternative source of genomic material for detection of specific mutations in certain types of cancers. EntroGen offers three ctDNA Mutation Detection Kits that allow non-invasive, ultra-sensitive testing to accurately detect the presence of mutations in ctDNA extracted from plasma.

ctDNA EGFR MUTATION DETECTION KIT

Targets	T790M, Ex19Del, L858R, C797S, L861Q, S768I, G719X, Ex20Ins	Sensitivity	As low as 10 copies
---------	--	-------------	---------------------



Low input requirements

ctDNA RAS MUTATION DETECTION KIT

Targets	KRAS G12, G13 NRAS G12, G13, Q61	Sensitivity	As low as 25 copies
---------	-------------------------------------	-------------	---------------------



Simple setup and interpretation

BRAF CODON 600 MUTATION ANALYSIS KIT II

Targets	BRAF V600 (E, E2, K, R, D, M, G)	Sensitivity	As low as 0.25%
---------	----------------------------------	-------------	-----------------



Compatible with most leading real-time PCR instruments

LEUKEMIA ASSAYS

Many hematologic malignancies carry characteristic chromosomal translocations believed to play an important role in leukemia pathogenesis. Detecting and identifying fusion gene transcripts in patient blood or bone marrow samples enables accurate diagnosis, prognosis, and treatment stratification for ALL, AML, CML and APL.

EntroGen's leukemia assays are all one-step RT-qPCR tests that offer results in under 2 hours with minimal hands-on time. Our assay portfolio includes a multiplexed translocation panel for initial diagnosis as well as separate kits targeting individual fusion gene transcripts for molecular monitoring during treatment.

Low positive controls are available for some leukemia assays.

2 QUANTIFICATION OPTIONS ABSOLUTE & RELATIVE		ALIGNED WITH THE INTERNATIONAL SCALE	
BCR-ABL P210 (Mbc _r) ONE-STEP DETECTION KIT			
Targets	b2a2, b3a2	Sensitivity	10 ^{-4.5} (0.0032%) For MRD

LEUKEMIA TRANSLOCATION PANEL			
Targets	AML1-ETO, CBFb-MYH11, MLL-AF4, TEL-AML1, E2A-PBX1, PML-RARA	Sensitivity	10 ⁻³ (0.1%) For MRD

AML1-ETO One-Step Detection Kit			
Targets	AML1-ETO	Sensitivity	10 ⁻³ (0.1%) For MRD

PML-RARA BCR1,BCR2,BCR3 One-Step Detection Kit			
Targets	PML-RARA	Sensitivity	10 ⁻⁴ (0.1%) For MRD

E2A-PBX1 One-Step Detection Kit			
Targets	E2A-PBX1	Sensitivity	10 ⁻³ (0.1%) For MRD

TEL-AML1 One-Step Detection Kit			
Targets	TEL-AML1	Sensitivity	10 ⁻³ (0.1%) For MRD

MLL-AF4 One-Step Detection Kit			
Targets	MLL-AF4	Sensitivity	10 ⁻³ (0.1%) For MRD

CBFB-MYH11 One-Step Detection Kit			
Targets	CBFB-MYH11	Sensitivity	10 ⁻³ (0.1%) For MRD

NPM1 ONE STEP MRD KIT			
Targets	Types A,B,D	Sensitivity	10 ^{-4.5} (0.0032%) For MRD

BCR-ABL P190 (mbcr) ONE-STEP DETECTION KIT			
Targets	e1a2 fusion transcript	Sensitivity	10 ⁻⁴ (0.01%) For MRD


IDH1/2 MUTATION DETECTION KIT			
Targets	IDH1 R100, R132 IDH2 R140, R172	Sensitivity	As low as 1% For MRD




NEXT GENERATION SEQUENCING

Next Generation Sequencing (NGS) is revolutionizing molecular diagnostics by advancing the capabilities of DNA sequencing far beyond traditional technologies. The speed, ultra-high throughput and multiplexing capabilities of NGS makes it possible to batch multiple samples, and interrogate hundreds of variants, all on a single run.

FEATURES

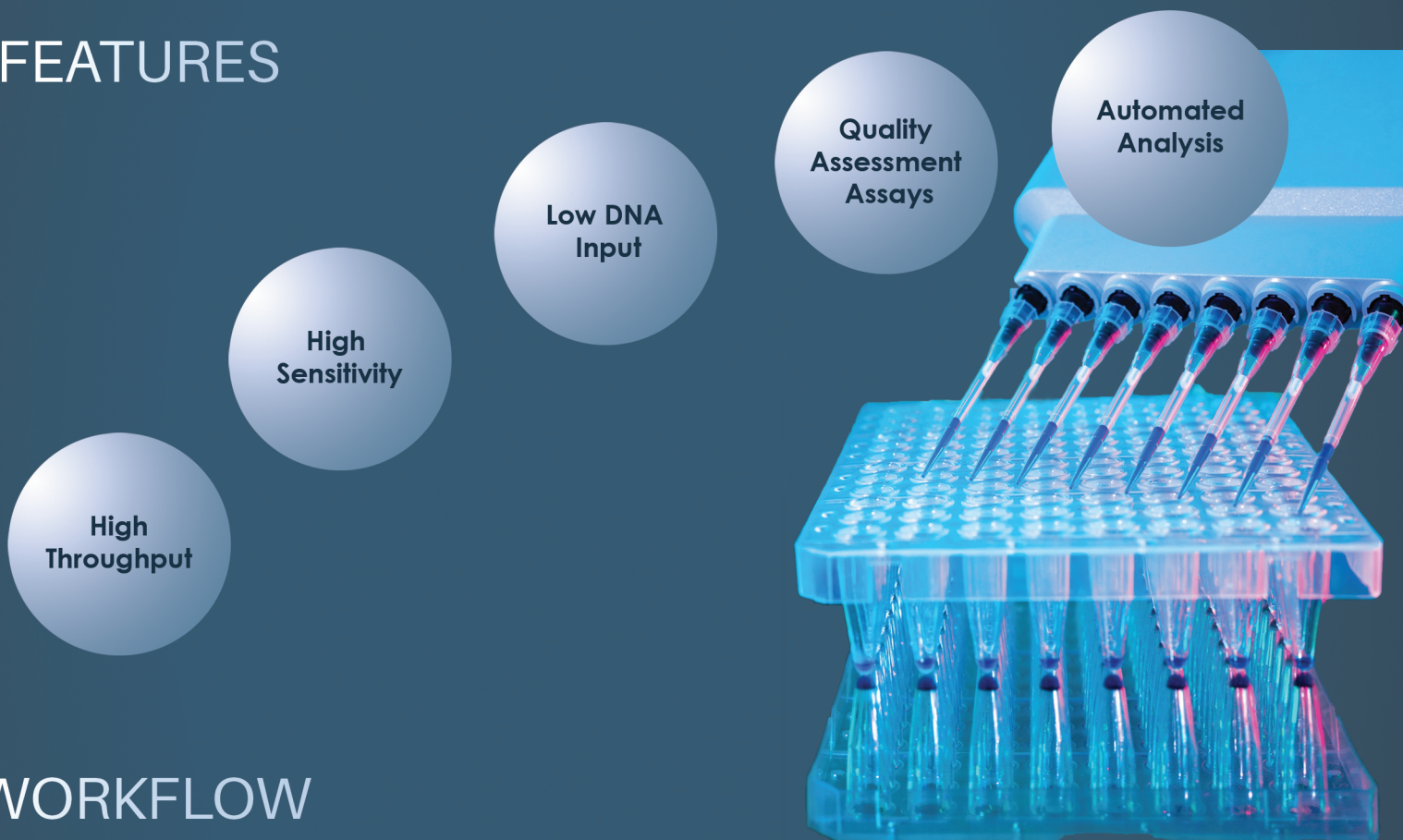

WORKFLOW

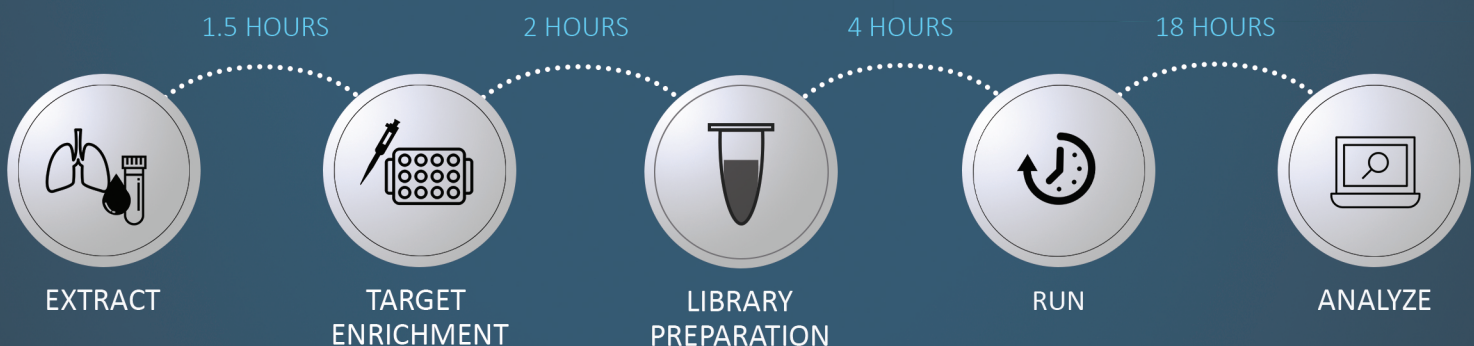

NGS ASSAYS

EntroGen NGS assays focus on identifying clinically relevant and actionable targets. A tiled amplicon PCR approach coupled with a high-fidelity polymerase for target enrichment are used to reduce allele dropouts and off target reads. Reagents for library preparation and clean-up are included, as well as user-friendly software that queries two independent databases for local data analysis.

BRCA COMPLETE EXPANDED PANEL FOR ILLUMINA® PLATFORMS

Targets	BRCA1/2, CHEK2, PALB2, RAD51C, TP53	Sensitivity	5%
Targets	POLE	Sensitivity	5%

Covers
Difficult-to-Sequence
TERT Promoter

NGS TARGETED HOTSPOT PANEL FOR ILLUMINA® PLATFORMS

Targets				Sensitivity	5%
BRAF	EXONS 11 & 15	KRAS	EXONS 2, 3, 4		
EGFR	EXONS 3, 7, 12, 15, 18-21	MET	EXONS 2, 11, 14, 16, 8-20		
ERBB2	EXONS 8, 14, 17-21, 24, 26	NRAS	EXONS 2, 3, 4		
HRAS	EXONS 2, 3	PDGFRA	EXONS 12, 14, 15, 18		
c-KIT	EXONS 2, 8-11, 13-15, 17, 18	PIK3CA	EXONS 1, 4, 6, 7, 9, 13, 18, 20		
IDH1	EXON 2	RET	EXON 10, 11, 13, 15, 16		
IDH2	EXON 4	TERT	PROMOTER		
JAK2	EXONS 12,14	TP53	ENTIRE EXOM		
POLE	EXONS 9-14	PTEN	EXONS 1-9	Sensitivity	5%

Detects
307 fusions in 23 genes

Boosted
NTRK Coverage

RNA FUSION GENE PANEL FOR ILLUMINA® PLATFORMS

Gene # of Detected Fusions								Sensitivity				5%			
ABL1	1	AXL	1	ERBB2	1	ETV4	1	FGFR2	10	NTRK1	34	PDGFRA	1	RET	41
ALK	61	BRAF	15	ERG	16	ETV5	3	FGFR3	26	NTRK2	18	PPARG	4	ROS1	35
AKT3	1	EGFR	2	ETV1	4	FGFR1	4	MET	10	NTRK3	16	RAF1	2		

SAMPLE QUALITY ASSESSMENT

INTERNAL QUALITY CONTROL ASSAY

Predicts optimal DNA input for maximal and reliable data output when using EntroGen's qPCR assays.

DNA FRAGMENTATION QUANTIFICATION ASSAY

Assesses sample quality by measuring the fragmentation ratio and quantifies the amount of amplifiable DNA for downstream qPCR and NGS assays.

LIBRARY QUANTIFICATION KIT FOR ILLUMINA®

Qualifies and quantifies library molecules that can be sequenced in downstream Illumina® assays.



Simple set-up and interpretation



Automated analysis available



Saves samples, reagents, and time by eliminating the need for reruns

MANUAL ISOLATION

pureNA® Cell-Free DNA Isolation Kit

DNA Fragment Size	100 - 500 bp	Reactions per Kit	50
-------------------	--------------	-------------------	----

pureNA® Genomic DNA Isolation Kit

Optimized for the removal of paraffin, partial reversal of formalin cross-linking and release of DNA from fixed tissue.	Reactions per Kit	50
---	-------------------	----



Magnetic bead-based technology



Simple, efficient workflow for DNA isolation



Compatible for downstream qPCR and NGS applications

DATA ANALYSIS SOFTWARE

EVA

EntroGen Variant Analysis Software

- Simplifies data interpretation and quickly determines:
 - variant calls
 - clinical relevance
 - current clinical trials
- Queries two independent databases for the most current genome annotations
- Complimentary, locally installed and available with all NGS assays



Import VCF File Sample ID: HD789A_S1 Search Mode: Standard Fast Annotate

View All Samples Clear Export Results Generate Report

HD789A_S1.vcf_085747335.vcf

Filters Clear [Export Current Page] [Go to page] Please select

Gene	Chr	Coord...	Type	Zyg...	Ref	Alt	Qual...	Filt...	DP	AD (Ref,Alt)	VF	Trans...	Cons...	HGVSc	HGVSp	dbSN...	ClinV...	ClinV...
BRCA2	chr13	32890572	SNV	Hete...	G	A	100	PASS	1013	788,225	0.2...	NM_0...	5 pri...	NM_000...		rs179...	VCV...	not
BRCA2	chr13	32890650	SNV	Hete...	G	A	100	PASS	1014	891,123	0.1...	NM_0...	misse...	NM_000...	NP_000...	rs803...	VCV...	beni
BRCA2	chr13	32893197	Inse...	Hete...	A	AT	72	R3x6	1624	1582,42	0.0...	NM_0...	intron...	NM_000...		rs276...	VCV...	unck
BRCA2	chr13	32893197	Dele...	Hete...	AT	A	100	R3x6	1631	1448,183	0.1...	NM_0...	intron...	NM_000...		rs276...	VCV...	unck
BRCA2	chr13	32905219	Inse...	Hete...	A	AT	100	R3x6	5896	5089,807	0.1...	NM_0...	intron...	NM_000...		rs751...	VCV...	beni
BRCA2	chr13	32905219	Inse...	Hete...	A	ATT	42	R3x...	5896	5805,91	0.0...	NM_0...	intron...	NM_000...		rs751...	VCV...	beni
BRCA2	chr13	32905219	Dele...	Hete...	AT	A	100	R3x6	5881	4762,1119	0.1...	NM_0...	intron...	NM_000...		rs751...	VCV...	beni
BRCA2	chr13	32905219	Dele...	Hete...	ATT	A	66	R3x...	5896	5794,102	0.0...	NM_0...	intron...	NM_000...		rs751...	VCV...	beni

Page 1 of 1, Found: 83, Total: 83

DNA Variant NGS Report

Variant Report


NT Change	AA Change	Variant Frequency	Zygosity	Type	Classification
BRCA2 (NM_000059.4)					
c.1114A>C	p.Asn372His	68.10%	Heterozygote	SNV	benign
c.4563A>G	p.Leu1521=	100.00%	Homozygote	SNV	benign
c.6513G>C	p.Val2171=	100.00%	Homozygote	SNV	benign
c.7397T>C	p.Val466Ala	99.90%	Homozygote	SNV	benign
c.7806-14T>C		50.40%	Heterozygote	SNV	benign
c.9257-16T>C		65.60%	Heterozygote	SNV	benign
c.9976A>T	p.Lys3326Ter	71.10%	Heterozygote	SNV	benign
TP53 (NM_000546.6)					
c.672+62A>G		99.90%	Homozygote	SNV	benign
c.376-91G>A		99.90%	Homozygote	SNV	benign
c.376-160_376-158del		67.10%	Heterozygote	deletion	benign, likely benign
c.376-283T>C		99.80%	Homozygote	SNV	benign
c.215C>G	p.Pro72Arg	99.60%	Heterozygote	SNV	benign
c.74-38C>G		99.90%	Homozygote	SNV	benign
BRCA1 (NM_007294.4)					
c.4535G>T	p.Ser1512Ile	63.50%	Heterozygote	SNV	benign
c.1067A>G	p.Gln356Arg	64.20%	Heterozygote	SNV	benign
c.442-34C>T		99.20%	Homozygote	SNV	benign
CHEK2 (NM_007194.4)					
c.1499C>G	p.Ser500Cys	30.50%	Heterozygote	SNV	uncertain significance
c.1149C>T	p.The383=	28.20%	Heterozygote	SNV	likely benign
c.1141A>G	p.Met381Val	23.50%	Heterozygote	SNV	uncertain significance

1. Variants without a ClinVar classification are not reported.
2. Filters applied:

DATA ANALYSIS SOFTWARE

EPAS

EntroGen PCR Analysis Software

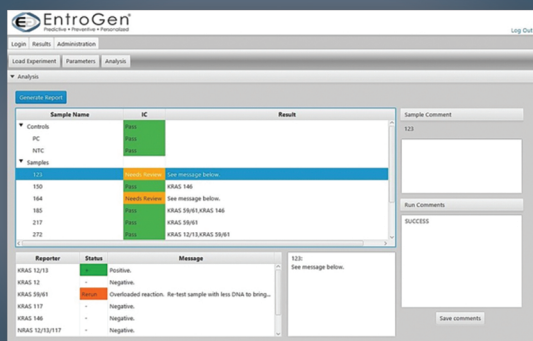

- Streamlines reporting and simplifies data interpretation
- Raw data files undergo QC check and mutation status call
- Simple point-and-click parameter setup and report generation
- Easy-to-interpret reports provide results that can be exported for integration with laboratory information management systems (LIS/LIMS)
- Complimentary, locally installed and available with select PCR assays

EntroGen Automation • Reagents • Software		Report Date: 10/16/2023 17:16			
Run Information:					
Run File Name:	BRC210_Run12345678.eds	Assay Name:	BCR210-ABS-QS5		
Lot No:	33835436 (Expiration: 08/2024)	Run Date:	02/27/2023 17:51:46.853		
Number of Samples:	4	User Name:	User Name		
EPAS / Configuration Version:	1.04 / 1.0	Instrument ID:	Instrument1234		
QC Info:					
Status:	Pass	Comments:			
Positive Control:	Pass				
Negative Control / NTC:	Pass				
Comments:					
Samples:					
Sample	IC Status	Results	Baseline Ratio	Log Reduction	Comments
MRD008	Pass	0.008	100	4.097	
MRD015	Pass	0.099	100	3.004	
MRD101	Pass	0.806	100	2.094	
MRD121	Pass	8.397	100	1.076	

BENCHbot⁴⁸⁰⁰

AUTOMATED SAMPLE PREP & PLATE SET-UP

Benchbot[®] 4800 Liquid Handling System is a fully automated sample preparation instrument that uses non-hazardous, eco-friendly reagents for nucleic acid isolation and PCR plate setup. Benchbot is designed to improve productivity and minimize hands-on time by prepping up to 48 samples simultaneously with the push of a button. The intuitive, user-friendly interface is simple to operate and comes pre-programmed with protocols for sample prep and PCR plate setup for EntroGen assays.

Features & Benefits



Simple, rapid, efficient workflow
from DNA/RNA isolation to plate setup



High throughput
Prep up to 48 samples simultaneously



Reduces hands-on time
No tedious manual reagent prep



Intuitive, user-friendly interface
Simple to operate - no advanced training required



Reduces risk of contamination




Reduces risk of human error



Delivers consistent, reproducible, reliable results



Easily integrated into any lab

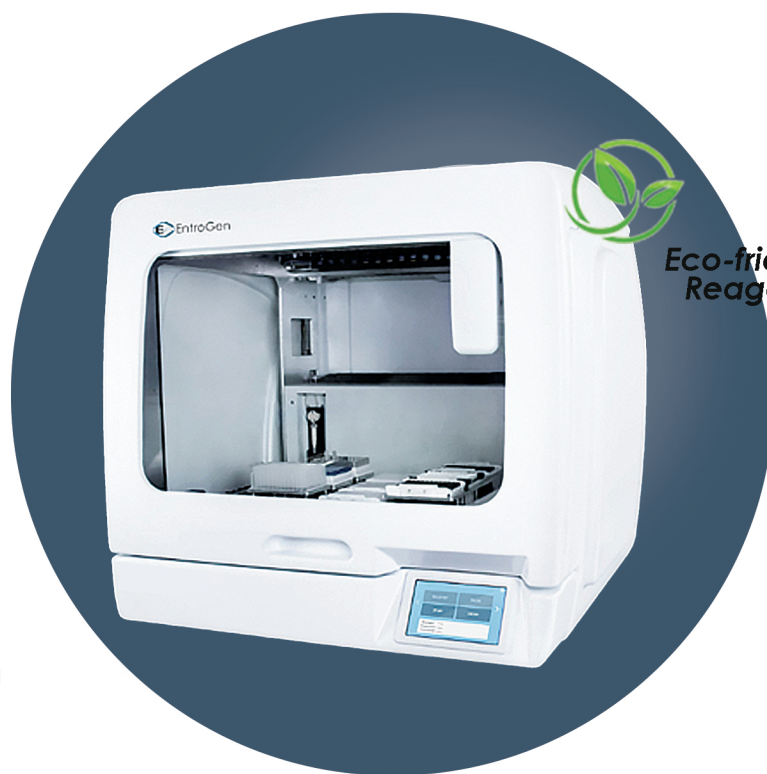

REAL-TIME PCR ASSAYS

COLORECTAL CANCER

Product Name	Catalog Number	Registration	RXN/Kit
KRAS MUTATION ANALYSIS KIT	KRAS-RT50	RUO, CE-IVD	50
KRAS/BRAF MUTATION ANALYSIS KIT	KRBR-RT50	RUO, CE-IVD	50
RAS MUTATION ANALYSIS KIT	RAS-RT50	RUO, CE-IVD	50
PIK3CA MUTATION ANALYSIS KIT	PI3K-RT48	RUO, CE-IVD	48
NRAS MUTATION ANALYSIS KIT	NRAS-RT50	RUO, CE-IVD	50
COLORECTAL CANCER MUTATION DETECTION PANEL	CRC-RT48	RUO, CE-IVD	48
CTDNA RAS MUTATION DETECTION KIT	CTRAS-RT40	RUO, CE-IVD	40
MLH1 METHYLATION DETECTION KIT	MLH1-RT24	RUO, CE-IVD	50
NTRK FUSION GENE DETECTION KIT	NTRK-RT48	RUO, CE-IVD	48
RASX MUTATION ANALYSIS EXTENSION KIT	RASX-RT50	RUO, CE-IVD	50
MICROSATELLITE INSTABILITY DETECTION KIT	MSI-RT24	RUO	24
CRCDX® RAS MUTATION DETECTION KIT	CRCDX-RAS	IVD	40

GASTROINTESTINAL STROMAL TUMORS

Product Name	Catalog Number	Registration	RXN/Kit
GIST MUTATION DETECTION KIT	GIST-RT44	RUO, CE-IVD	44

GLIOBLASTOMA

Product Name	Catalog Number	Registration	RXN/Kit
MGMT METHYLATION DETECTION KIT	MGMTX-RT44	RUO, CE-IVD	60
IDH1/2 MUTATION DETECTION KIT	IDH-RT38	RUO, CE-IVD	38

LUNG CANCER

Product Name	Catalog Number	Registration	RXN/Kit
CTDNA EGFR MUTATION DETECTION KIT	CTEGFR-RT48	RUO, CE-IVD	48
EGFR MUTATION ANALYSIS KIT	EGFR-RT52	RUO, CE-IVD	52
EGFRX MUTATION ANALYSIS EXTENSION KIT	EGFRX-RT52	RUO, CE-IVD	52
LUNG CANCER RNA PANEL	LUNG-RT48	RUO, CE-IVD	48

MELANOMA

Product Name	Catalog Number	Registration	RXN/Kit
BRAF CODON 600 MUTATION ANALYSIS KIT	BRAFX-RT64	RUO, CE-IVD	64
C-KIT MUTATION DETECTION KIT	CKIT-RT44	RUO, CE-IVD	44

THYROID CANCER

Product Name	Catalog Number	Registration	RXN/Kit
THYROID CANCER MUTATION KIT	THDNA-RT64	RUO, CE-IVD	64
THYROID CANCER FUSION GENE DETECTION KIT	THRNA-RT32	RUO, CE-IVD	32

LEUKEMIA MRD ASSAYS			
Product Name	Catalog Number	Registration	RXN/Kit
BCR-ABL P190 (mbcr) ONE STEP DETECTION KIT	BCR190-QRT46	RUO, CE-IVD	46
BCR-ABL P210 (Mbcr) ONE STEP DETECTION KIT	BCR210-QRT46	RUO, CE-IVD	46
LEUKEMIA TRANSLOCATION PANEL V2	LEUKMP-RT24	RUO, CE-IVD	24
AML1-ETO ONE STEP DETECTION KIT	LEUK1-QRT24	RUO, CE-IVD	24
E2A-PBX1 ONE STEP DETECTION KIT	LEUK2-QRT24	RUO, CE-IVD	24
MLL-AF4 ONE STEP DETECTION KIT	LEUK3-QRT24	RUO, CE-IVD	24
PML-RARA bcr1, 2, 3 ONE STEP DETECTION KIT	LEUK4-QRT24	RUO, CE-IVD	24
TEL-AML1 ONE STEP DETECTION KIT	LEUK5-QRT24	RUO, CE-IVD	24
CBFb-MYH11 ONE STEP DETECTION KIT	LEUK6-QRT24	RUO, CE-IVD	24
NPM1 ONE STEP MRD KIT	NPM1-QRT24	RUO	24

NGS ASSAYS

BREAST CANCER			
Product Name	Catalog Number	Registration	RXN/Kit
BRCA COMPLETE™ EXPANDED PANEL	BRCAX-NGS	RUO, CE-IVD	40

FUSION PANEL			
Product Name	Catalog Number	Registration	RXN/Kit
RNA FUSION GENE PANEL	RFGP-NGS	RUO, CE-IVD	40

HOTSPOT PANEL			
Product Name	Catalog Number	Registration	RXN/Kit
NGS TARGETED HOTSPOT PANEL	THSP-ILNGS	RUO, CE-IVD	40

NGS ADDONS			
Product Name	Catalog Number	Registration	RXN/Kit
POLE HOTSPOT NGS ADDON	NGS-POLE	RUO	40
PTEN HOTSPOT NGS ADDON	NGS-PTEN	RUO	40

ISOLATION KITS, SAMPLE QUALITY ASSESSMENT & REFERENCE CONTROLS

Product Name	Catalog Number	Registration	RXN/Kit
pureNA® CELL FREE DNA ISOLATION KIT	CFDNA-50	RUO, CE-IVD	50
pureNA® GENOMIC DNA ISOLATION KIT	GDNA-50	RUO, CE-IVD	50
INTERNAL QUALITY CONTROL ASSAY	IQCA-RT50	RUO, CE-IVD	50
DNA FRAGMENTATION QUANTIFICATION ASSAY	FQA-RT40	RUO, CE-IVD	40
LIBRARY QUANTIFICATION KIT FOR ILLUMINA®	LIBQ-NGS	RUO, CE-IVD	75
5% MGMT METHYLATION REFERENCE CONTROL	MGMT-LPC	RUO	N/A
LOW POSITIVE CONTROL FOR LEUKMP-RT24	LEUKMP-LPC	RUO	N/A
LOW POSITIVE CONTROL FOR LEUK4-RT24	LEUK4-LPC	RUO	N/A
LOW POSITIVE CONTROL FOR BCR190-QRT46	BCR190-LPC	RUO	N/A
LOW POSITIVE CONTROL FOR BCR210-QRT46	BCR210-LPC	RUO	N/A
BCR-ABL P210 QUANTITATIVE CONTROL STANDARDS	BCR210-AQSTD	RUO	N/A





Automation



Reagents



Software

WWW.ENTROGEN.COM



INFORMATION@ENTROGEN.COM

MB-10142024

© Copyright 2024, EntroGen Inc. All rights reserved.

EntroGen products are for Research Use
Only (RUO) in the US and CE-IVD in Europe.